

MedibinoNEO[®]: Baby positioning aid for clinical use



medibino^{neo}



KLUBA
MEDICAL

Positioning aid & pressure relief

Give every baby the benefit of gentle but stable positioning

Safe, developmentally supportive and durable head positioning system and stabilisation aid for prematures and newborns

Soft, comfortable support for the baby's head and reliable positioning aid for mid-head position and therapeutic interventions

- Helps minimise head deformation and skin abrasion
- Assists in stabilising the midline position of the head
- Helps relieve pressure and promotes correct positioning
- Provides support for the soft skull of the premature baby
- The tested gel is guaranteed free of mould and bacteria
- Easy cleaning of the gel ring - use for up to one year under normal conditions
- Different disposable and multi-weave covers guarantee flexible use according to the needs of the premature baby unit
- All materials are tested and selected according to the soft skin of premature babies

3 Sizes Meeting all developmental needs



MedibinoNEO® size S for premature babies up to 800 g,
Dimensions:
90 mm x 80 mm x 18 mm


MedibinoNEO® size M for premature babies between approx. 800g and 2000g,
Dimensions:
110 mm x 95 mm x 22 mm

MedibinoNEO® size L for premature babies over 2000g and newborns,
Dimensions:
130 mm x 110 mm x 25 mm

Disposable, baby-soft covers available for each size



Disposable soft cover

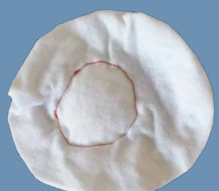

Disposable soft cover WITH gel ball depot included



Soft but stable gel ring helps avoid frequent repositioning in supine or side position

Three sizes allow flexible use for all babies


Inner cushion filled with gel balls provides additional stability when needed

Fresh colours guarantee retrieval and provide a child-friendly design

With MedibinoNEO[®], every baby can benefit from comfortable and long-lasting positioning support and pressure relief

The gel ring is an effective way to minimise pressure and protect the soft skull of the premature baby. The ring shape minimises the need for repositioning in the supine and lateral positions and helps relieve pressure points. For additional stabilisation, the MedibinoNEO[®] has a soft inner cushion with gel balls. The soft surface of the different covers makes babies feel more comfortable and allows flexible use depending on the needs of the clinic.

MedibinoNEO[®] is an effective helper for various applications on the ward

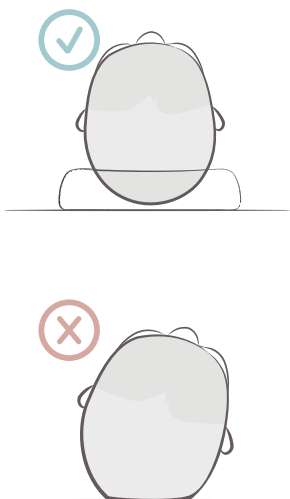

Apply MedibinoNEO for:

- For the pressure-relieving positioning of premature babies and infants during a long stay in hospital as well as in the maternity ward to inform parents about measures to prevent cranial deformities
- For positioning in 30° elevation and for stabilisation of the central head position in asphyxia and, if required, for prevention of intraventricular haemorrhages (cerebral haemorrhage)
- To stabilise the head in the central position during necessary interventions

Product Specifications

1. Soft gel storage ring made of TG Pur Gel; outer material: biocompatible PU film

The inner core of MedibinoNEO[®] is non-sterile. It can be disinfected by wiping and is suitable for cleaning with commercially available disinfectants. Shelf life is 24 months; application period from first use is 12 months.

2. Inner cushion coated, filled with soft granules

The inner pad is non-sterile. It is water-repellent and can be disinfected by wiping and is suitable for cleaning with commercially available disinfectants. The shelf life is 24 months; the application period from the first use is 12 months.

3. Cover made of GOTS cotton, medically tested

The cover of the storage ring can be washed at 60°C.

4. Disposable covers made of cotton knit, medically tested

The disposable covers are available WITH and WITHOUT soft granules and they are designed for single use.

Sizes & Components

MedibinoNEO ring – TG PUR Gel, outer material biocompatible PU film

24021703-L	size L, 130 mm x 110 mm x 25 mm	6 pcs. per box
24021703-M	size M, 110 mm x 95 mm x 22 mm	6 pcs. per box
24021703-S	size S, 90 mm x 80 mm x 18 mm	6 pcs. per box

MedibinoNEO inner pillow, WWL 200 filled with soft granulat (optional)

25021703-L	size L, for use with ring size L	6 pcs. per box
25021703-M	size M, for use with ring size M	6 pcs. per box
25021703-S	size S, for use with ring size L	6 pcs. per box

MedibinoNEO cover GOTS organic cotton with hotel closure

26021703-L	size L, for use with support ring and inner cushion (optional) size L	6 pcs. per box
26021703-M	size M, for use with support ring and inner cushion (optional) size M	6 pcs. per box
26021703-S	size S, for use with support ring and inner cushion (optional) size S	6 pcs. per box

MedibinoNEO disposable Cover cotton knit WITHOUT soft granules - colour coded

27021703-L	size L, for use with support ring and inner cushion (optional) size L	50 pcs. in polybag
27021703-M	size M, for use with support ring and inner cushion (optional) size M	50 pcs. in polybag
27021703-S	size S, for use with support ring and inner cushion (optional) size S	50 pcs. in polybag


MedibinoNEO disposable cover cotton knit WITH soft granules - colour coded

31021703-L	size L, for use with ring size L	50 pcs. in polybag
31021703-M	size M, for use with ring size M	50 pcs. in polybag
31021703-S	size S, for use with ring size L	50 pcs. in polybag




Kluba Medical GmbH
Brunnenstraße 23
40223 Düsseldorf
+49 (0) 211 5421 2590
info@kluba-medical.com

Medical Background

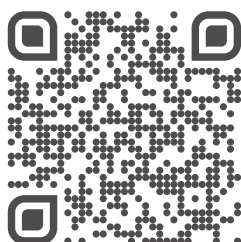

Click here or Scan

Instructions for use

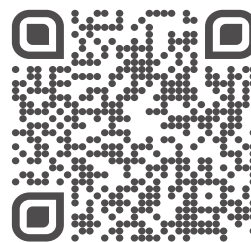

Click here of Scan

Fields of application



Click here of Scan